

YA-M5-6

Manual hospital bed



*Better care
Better comfort!*

YA-M5-6

MANUAL HOSPITAL BED

DESCRIPTION

The YA-M5-6 is a 5 functions manual hospital bed with ergonomic design elegant appearance. Patient safety and functionality at low cost are the priority concerns for this bed. It provides a complete range of bed postures with low budget.

- FEATURES
- Bed Frame electrophoresis and powder coating
 - MS perforated platform with ventilation hole, molded
 - Four split side rail, drop down under the bed platform
 - Head and Footboard Detachable head/footboard with easy lock
 - 4 pieces foldable hand cranks
 - urine bag hooks on both sides
 - With angle indicator embeded in the railing
 - Oxygen bottle holder and IV pole

STANDARD SPECIFICATION

Overall size	2240*1085*450-700mm
Recommended mattress size	2015*875mm
Height adjustable	450-700mm
Backrest adjustable	0° to 70°
Footrest adjustable	0° to 40°
Trendelenburg/reverse trendelenburg	0-12°
Castor	Dia 125mm
Loading capacity	200kg



Section adjustable

Back, foot section and height can be adjusted by crank separately, can meet the needs of patients in different positions.



Trendelenburg/reverse trendelenburg

Trendelenburg/reverse trendelenburg are adjusted by manual crank.



Manual crank

Four pieces foldable manual crank at foot side of the bed, and there is a medical record holder at the end of the bed for easy identification by medical staff.



Angle display

Back, foot section and height can be adjusted by crank separately, can meet the needs of patients in different positions.




IV hole

There are infusion holes on the four corners of the bed.



Urine bag hook

Urine bag hooks on both sides of the bed for easy use in any position.

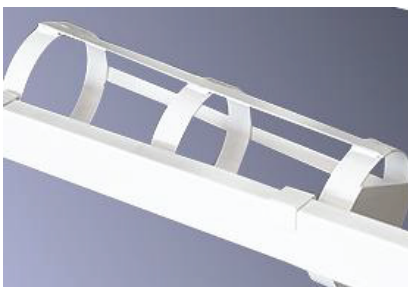

Oxygen bottle holder

There is an oxygen bottle holder at the bottom of the bed for easy use.



Side rail lock

Tuck-away side railing meets IEC60601-2-52 safety standard, special protection mechanism can prevent patient to drop down the railing from inside bed.



Castors

4 luxurious silent double-sided wheels with individual brake.